
waterpik®

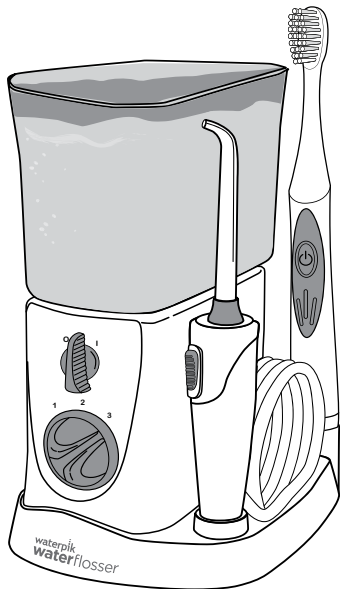Waterpik® Water Flosser plus Sonic Toothbrush Model WP-700

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www.waterpik.com

TABLE OF CONTENTS

Important Safeguards	3
Product Description	5
Getting Started With Your Water Flosser	7
Cleaning and Troubleshooting Your Waterpik® Water Flosser	10
Getting Started With Your Sonic Toothbrush/ Care of Your Sonic Toothbrush	12
Limited Two-Year Warranty	13

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. or your local distributor which can be found at www.waterpik.com/intl.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips, brush heads, and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

- Do not use iodine, bleach, or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- If the Pik Pocket® tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket® tip.
- Do not use this product for more than 5 minutes in each two-hour period.
- Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.
- There are no consumer serviceable parts in this appliance and it does not require maintenance.
- Fill reservoir with warm water or other dental professional recommended solution only.


Battery Notice (Toothbrush)

- This toothbrush runs on 1 AAA alkaline battery. For optimum performance, replace battery with the same alkaline type.
- Remove exhausted battery from the toothbrush and dispose of properly.

NOTE: Never leave an exhausted battery in your toothbrush.

SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

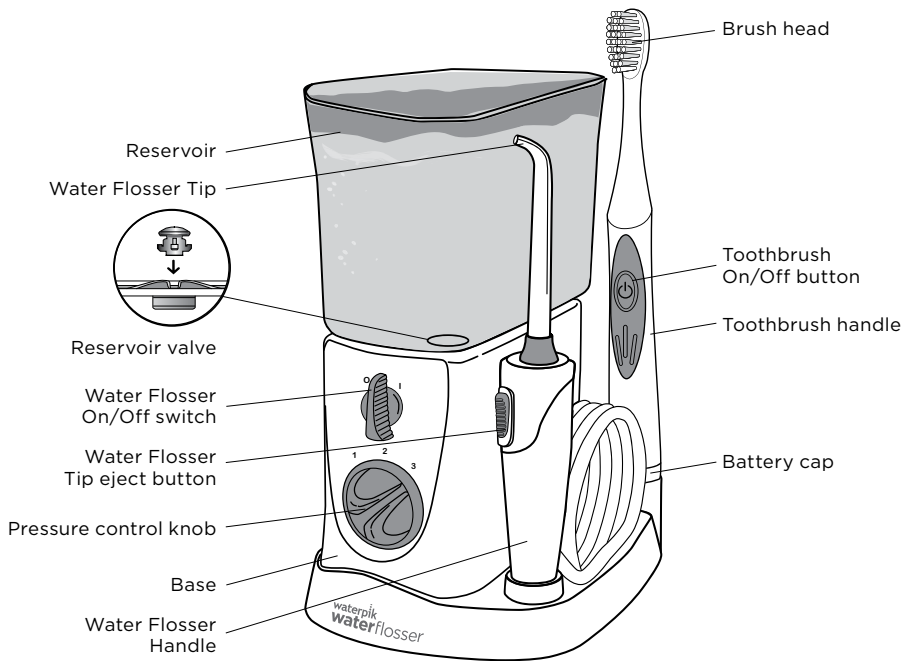Replace Every 3 Months	
Orthodontic Tip	BEST FOR
	☒ Braces ☒ General Use
Pik Pocket® Tip	BEST FOR
	☒ Periodontal Pockets ☒ Furcations
Plaque Seeker® Tip	BEST FOR
	☒ Implants ☒ Bridges ☒ Crowns ☒ General Use
Standard Brush Head	BEST FOR
	☒ Plaque Removal

Replace Every 6 Months	
Classic Jet Tip	BEST FOR
	☒ General Use
Tongue Cleaner	BEST FOR
	☒ Fresher Breath


NOTE: All tips not included in all models.

Replacement tips/attachments may be purchased from your local distributor, which can be found at www.waterpik.com/intl.


GETTING STARTED WITH YOUR WATER FLOSSER

Inserting and Removing Tips

Insert tip into the center of the Water Flosser handle. The colored ring will be flush with the top of the handle if the tip is correctly locked into place.

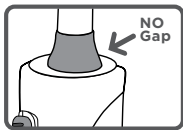To remove tip from the handle slide the tip eject button up towards the top of the handle and lift tip from handle.

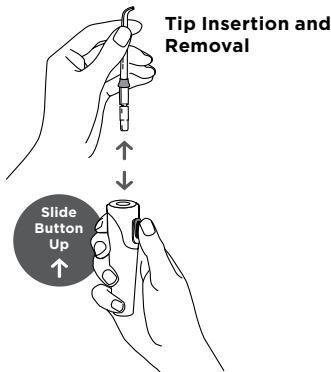
Powering the Unit

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.


Preparing the Reservoir

Remove the reservoir and fill with lukewarm water. Replace the reservoir over the base and **press down firmly**.

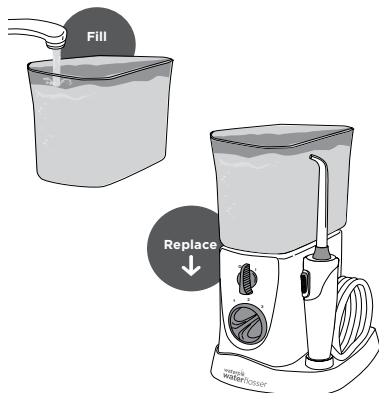
Prime Unit (First Time Use Only)

Set pressure control dial on base to high, point the handle and tip into sink. Turn unit ON (I) until water flows from the tip. Turn unit OFF (O).

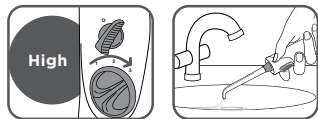Adjusting the Pressure Setting

Turn the pressure control dial on the base of the unit to the lowest setting (pressure setting No. 1) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

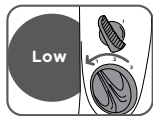Recommended Technique

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.


Lean over sink. Direct the stream at a 90-degree angle to your gumline with tip in mouth.

Slightly close lips to avoid splashing but allow water to flow freely from mouth. Set pressure at desired setting and turn the unit ON (I). Begin cleaning.

When Finished

Turn the unit OFF (O). Empty any liquid left in the reservoir. Unplug the unit.

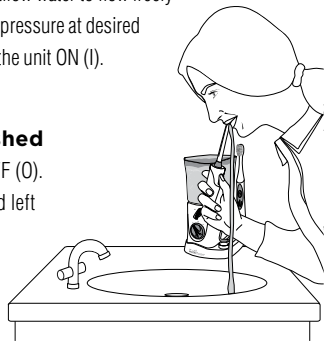
TIP USAGE

Pik Pocket® Tip


The Pik Pocket® Tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket® Tip, **set the unit to the lowest pressure setting.** Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

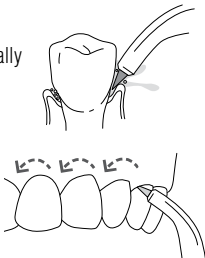
Tongue Cleaner Tip

To use the Tongue Cleaner Tip, set to lowest pressure setting, place tip in the center/middle of your tongue about half way back and start the Water Flosser. Pull forward with light pressure. Increase pressure as you prefer.

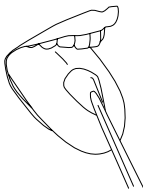Orthodontic Tip

To use the Orthodontic Tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

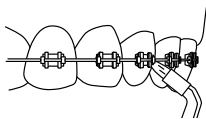Plaque Seeker® Tip


To use the Plaque Seeker® Tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

CLEANING AND TROUBLESHOOTING YOUR WATERPIK® WATER FLOSSER

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

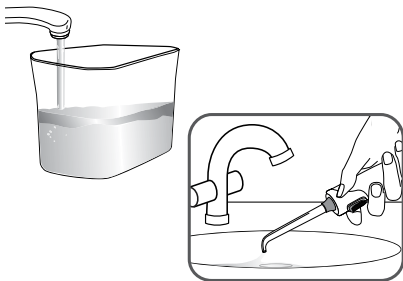Cleaning

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe.

NOTE: Remove reservoir valve before placing reservoir in dishwasher.

Remove

Replace

Removing Hard Water Deposits

Add 2 tablespoons of white vinegar to a full reservoir of warm water. Set pressure control dial to high, point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

Leakage

If leaking due to a damaged hose, replacement hoses and handles are available from your local distributor which can be found at www.waterpik.com/intl.

If leakage occurs from the reservoir, ensure the valve is right side up or contact your local distributor, which can be found at www.waterpik.com/intl, for a replacement valve.

Refer to the serial and model numbers in all correspondence. These numbers are located on the back of the handle and bottom of the reservoir respectively.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com/intl to find your local distributor. Please do not ship product to Waterpik International, Inc. offices, this will delay service.

Some replacement parts and accessories for the Water Flosser are also available from your local distributor which can be found at www.waterpik.com/intl.

GETTING STARTED WITH YOUR SONIC TOOTHBRUSH/ CARE OF YOUR SONIC TOOTHBRUSH


Battery insertion

To insert the battery:

1. Unscrew the battery cap.

2. Insert the AAA alkaline battery into brush with positive end first.

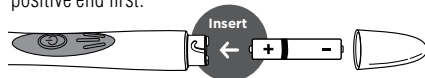3. Screw on the battery cap. Replace the AAA alkaline battery with a new AAA alkaline battery when it is exhausted.

Sonic brushing

- Slide the brush head onto the handle. Push the brush head down until it is fully seated on the handle.
- Place the brush head into your mouth before turning the unit on.
- For best results, position brush along gumline at a 45-degree angle toward the gum.

- Use light pressure with a gentle back and forth motion covering one or two teeth at a time. Overlap the teeth as you work your way around the mouth.
- Move the brush slowly around entire mouth.

Cleaning

1. After each use, turn off unit, and rinse as you would a manual toothbrush.
2. Wipe the handle dry after each use.

Service

Other than the replaceable battery, the sonic toothbrush has no consumer serviceable electrical items and does not require routine service maintenance. For all your service needs go to www.waterpik.com/intl to find your local distributor. Please do not ship product to Waterpik International, Inc. offices, this will delay service.

Battery Disposal

Please dispose of the replaceable battery per local disposal guidelines and/or laws.

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found on [www.waterpik.com/intl]. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence.

To locate your local distributor, go to www.waterpik.com/intl.

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Water Pik, Inc.

1730 East Prospect Road
Fort Collins, CO 80553-0001 USA
www.waterpik.com

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Brazil, Canada, Chile, China, Czech Republic, Denmark, EU, Finland, France, Germany, Hong Kong, Hungary, India, Italy, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Poland, Russian Federation, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Turkey, Ukraine, United Kingdom, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Japan, Mexico, Korea, Russian Federation, Switzerland, Trinidad & Tobago, and the United States.

Pik Pocket® is a trademark of Water Pik, Inc. registered in Canada and the United States.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

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Form No. 20016966-F AB
©2014 Water Pik, Inc.

Waterpik, International, Inc.

Castle Court, 41 London Road
Reigate, Surrey RH2 9RJ
United Kingdom

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Pik Pocket® è un marchio registrato di Water Pik Inc. registrato in Canada e negli Stati Uniti.
Plaque Seeker® è un marchio di Water Pik, Inc. registrato in Canada, Cina, UE e negli Stati Uniti.
Pik Pocket® ist eine Marke von Water Pik, Inc., registriert in Kanada und in den Vereinigten Staaten.
Pik Pocket® is een handelsmerk van Water Pik, Inc., geregistreerd in Canada en de Verenigde Staten.
Plaque Seeker® is een handelsmerk van Water Pik, Inc. geregistreerd in Canada, China, EU, en de Verenigde Staten.
Pik Pocket® est une marque de commerce de Water Pik, Inc. déposée au Canada et aux États-Unis.